



Osteo-x[®]

External fixation, distraction and reduction of small bones



- Simple application
- Compact and lightweight to ensure comfort for the patient
- Manufactured from PEEK and carbon fibre
- Radiolucent to aid viewing of the fracture site
- Pre-assembled and sterile, ready for immediate use

The Osteo-x[®] external fixator is a device used for the fixation and distraction of small bones, primarily in the hand. It incorporates four K-wires which are driven into the bone and placed to provide three-dimensional stability. The displacement between the blocks can be set prior to attachment and the distraction/compression adjusted by rotating the threaded shaft with the supplied tool.

The simple nature of the device reduces operative time as it is supplied fully assembled and sterile packed.

Product


Catalogue No.	Description
OSTX-ADJ-SP	Osteo-x [®] Sterile Adjustment Tool
OSTX-ASM-SP	Osteo-x [®] External Mini Fixator



Related Products

ChiroKlip[®] - Protecting the protruding ends of K-wires

Catalogue No.	Description
ChiroKlip [®]	ChiroKlip [®] - Sterile - 1 box (50 per box)
Quantity discounts available	

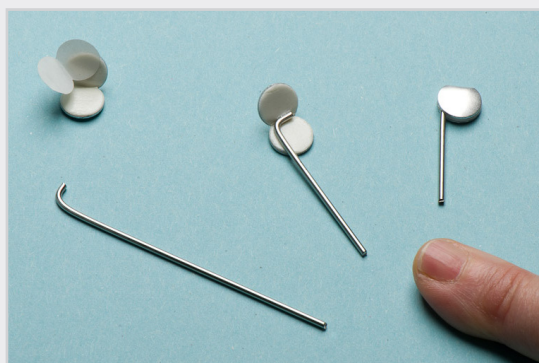

Osteotec Limited

9 Silver Business Park, Airfield Way, Christchurch, Dorset, BH23 3TA
 Tel: +44 (0)1202 487 885 Fax: +44 (0)1202 487 886
 Web: www.osteotec.co.uk Email: sales@osteotec.co.uk